# Quick Start Guide


## GETTING STARTED

- 1 With the unit turned OFF(O), connect the charging cable to user supplied USB wall charger. Plug the charger into the wall outlet, then connect the other end of the charging cable to the unit, and charge.

**Note:** Refer to the IMPORTANT SAFEGUARDS section for required USB wall charger specifications.


Make sure unit is dry before opening charge port cover and connecting the charger. Replace the charge port cover after charging and ensure it's fully seated (flush with the outside surface of the handle).

**Charge Unit for 4 Hours Prior to First Use.**

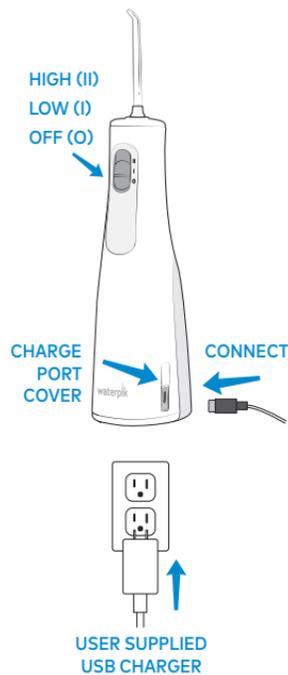
- 2 Press tip firmly into unit until it clicks into place. To remove tip, push the tip eject button.

**Caution: Do not press tip eject button while unit is running.**

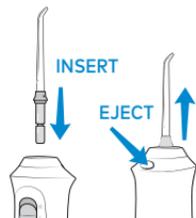- 3 Unplug USB charging cable from the unit. Ensure charge port cover is in place prior to filling with water and during use. Lift the flip top on reservoir and fill with lukewarm water.

Refill as needed.

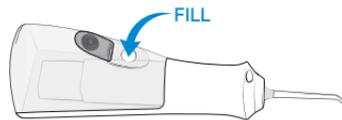
## WATER FLOSS

- 4 Lean over sink and close lips enough to prevent splashing while still allowing water to flow from mouth into the sink.

Keep unit upright during use to ensure continuous water flow.


- 5 For first time use, adjust the pressure setting to the LOW(L) setting.
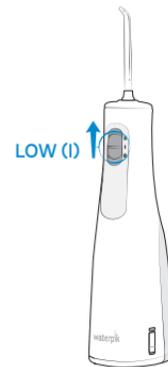
## TECHNIQUE

- 6 Aim water at the gumline at a 90 degree angle.

Follow the gumline and pause briefly between teeth.

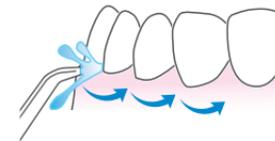**TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.**

## CONTACT US

To watch a How to Video and find additional product information, visit our website:

Canada: [waterpik.ca/enhance](https://www.waterpik.ca/enhance)  
Call Us! 1-888-226-3042 (CAN)

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