waterpik®

Congratulations!

You are now on your way to better oral health!

✓ Register and Save

Register your product at waterpik.ca/registration to receive over \$50 in savings (includes 50% off your next order of assorted tips + tip storage case).

✓ We're Here to Help

We hope you are loving your new Cordless Express Water Flosser. If for any reason you are not satisfied or have any product questions, please contact our fast and friendly Customer Service Team based in Fort Collins, Colorado.

The Cordless Express Water Flosser is backed by a 1-year manufacturer's warranty.

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Register and Save!

Save over \$50!

Register at Waterpik.ca/Registration to receive exclusive discount code.

Offer includes **50% off your next order** of assorted tips + tip storage case.

We'd Love to Hear from You!

Connect with Us
[@Waterpik](https://www.instagram.com/waterpik)

Waterpik® 14 Day Challenge App

Make water flossing a habit! Try our new mobile App for the best water flossing experience.

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
GETTING STARTED / POUR COMMENCER

- 1** Remove reservoir

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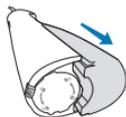- 2** To remove battery door, twist, and pull.

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- 3** Insert 3 AA batteries.

Insérez 3 piles AA.

- 4** To replace battery door, line-up arrows, push and twist.

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- 5** Insert tip firmly until it clicks into place.

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- 6** If the reservoir was removed, insert reservoir back on the unit. Fill reservoir with warm water.

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WATER FLOSS / HYDROPULSEUR

- 7** With water switch in OFF position, place tip in mouth. Slide power switch to low (I) and increase to your liking.

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- 8** Lean over sink. Close lips enough to prevent splashing. Keep unit upright during use to ensure continuous water flow.

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TECHNIQUE / TECHNIQUE

- 9** Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.

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NEED HELP? / BESOIN D'AIDE?

Visit Waterpik.ca to watch a How To Video and find additional product support.

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