
**waterpik**<sup>TM/MC</sup>
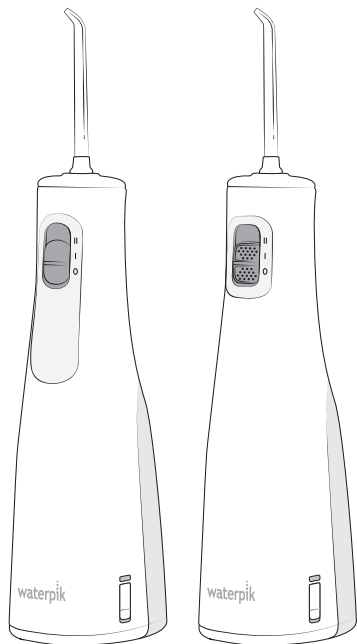
# **Waterpik<sup>TM</sup> Water Flosser WF-20 | 21 Series**

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[www.waterpik.ca](http://www.waterpik.ca)

# IMPORTANT SAFEGUARDS

## READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

## DANGER:

### To reduce the risk of electrocution:

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not handle the user supplied USB wall charger with wet hands.
- Do not place or drop the device or user supplied USB wall charger into water or other liquid.
- Check the charging cord and user supplied USB wall charger for damage before the first use and during the life of the product.

## WARNING:

### To reduce the risk of burns, electrocution, fire, or injury to persons:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories recommended by Water Pik, Inc.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc. at 1-888-226-3042 or visit [www.waterpik.ca](http://www.waterpik.ca).
- Keep the charging cord and user supplied USB wall charger away from heated surfaces.
- Never drop or insert any object into any opening or hose.

- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not expose the battery or appliance to excessive temperatures.
- A user supplied USB-A wall charger (power adapter), safety certified (IPX4, Class II) and with a rated output of 5.0 VDC, 1.0 A Max must be used to charge this device.
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Unplug USB wall charger when not in use.
- Do not use the Water Flosser for more than five minutes in each two hour period.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into the nose or ear.
- If your physician advises you to receive antibiotic premedication before dental procedures, you should consult your dental professional before using this product or any other oral hygiene aid.
- Fill reservoir with warm water or other dental professional recommended solutions only.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce performance and will shorten the life of the product.
- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with a new Pik Pocket™ tip.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- Do not use the Water Flosser for more than five minutes in each two hour period.
- There are no consumer serviceable electrical parts in this product and it does not require maintenance.
- Remove charging cord from the unit and replace charging port cover prior to filling the reservoir.
- If the supply cord is damaged, it must be replaced by a special cord or assembly available from the manufacturer or its service agent.

**SAVE THESE  
INSTRUCTIONS**

## TABLE OF CONTENTS


Product Description	5
Getting Started	6
Cleaning and Troubleshooting Your Waterpik™ Water Flosser	9
Battery Care, Battery Disposal/Recycling, Battery Removal	10
Limited Two Year Warranty	11

### SYMBOLS

The following symbols may appear on your product and/or packaging. Some of these symbols may not be relevant in your region and are listed for information purposes only.

Symbol	Definition	Symbol	Definition	Symbol	Definition	Symbol	Definition
	Refer to instruction manual		Detachable Power Supply		Do not use this device in a shower, bathtub and do not immerse in any liquid.		Regulatory Compliance Mark
	Manufacturer		Direct Current (DC) Voltage		Authorized Representative in the European Community		Separate collection for waste electrical and electronic equipment.
IPX7 or IPX4	Ingress Protection Code		Alternating Current (AC) Voltage		CE marking of conformity		
	IEC Class II Equipment		Do not immerse this device in any liquid.		UK Conformity Assessed		

1. Power handle
2. Power/pressure switch
3. Battery charge indicator
4. Charge port/cover
5. Tip
6. Tip eject button
7. Reservoir fill door
8. Reservoir
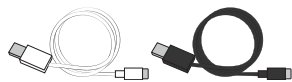
#### USB Charging Cable


20032330-2 White  
20032330-4 Black

WF-20

WF-21
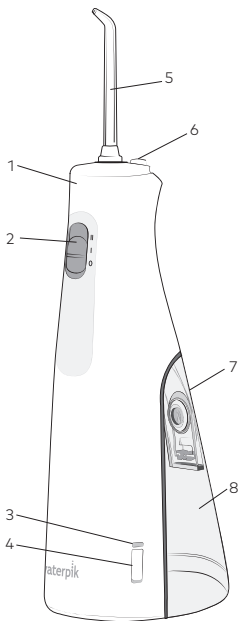
#### TIPS

##### A. Orthodontic Tip\*

- Braces/General use

##### B. Pik Pocket™ Tip\*

- Periodontal Pockets/Furcations

##### C. Plaque Seeker™ Tip\*

- Implants/Crowns/Bridges/  
Retainers/General use

##### D. Classic Jet Tip\*\*

- General use





##### E. Tongue Cleaner\*\*

- Fresher breath

##### F. Travel Plug

- Keeps excess water from  
draining out while traveling


#### NOT ALL TIPS INCLUDED IN ALL MODELS

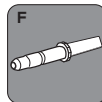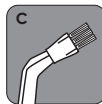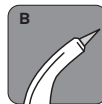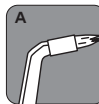\*Replace Every 3 months \*\*Replace Every 6 months

Replacement tips/attachments may be purchased online at  
[www.waterpik.ca](http://www.waterpik.ca), or by phone from Water Pik, Inc. at 1-888-226-3042.

## GETTING STARTED

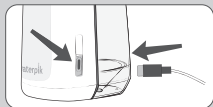USB wall charger is not standard accessory with this appliance.

### Charging Unit Prior to First Use

Connect the USB charging cable to the unit (USB-C) and then to the USB wall charger (USB-A).

Then plug the USB wall charger into the power outlet and charge unit for approximately 4 hours.

**Note: Refer to the IMPORTANT SAFEGUARDS section for required USB wall charger specifications.**

**Always ensure USB plugs and USB ports are clean and dry before connecting.**

LED will glow as it charges. LED stays lit when fully charged until USB cable is unplugged from the device. LED blinks 8 times once you are finished flossing and if the battery is below 30%.

**NOTE: Make sure unit is dry before opening weathercap and connecting the charger.**

**Replace the weathercap after charging and ensure it's fully seated (flush with the outside surface of the handle).**

### Filling the Reservoir

Before filling, disconnect charging cable and replace charging port cover.


Lift flip top on reservoir and fill with lukewarm water.


### Inserting and Removing Tips

Insert tip into the center of the knob at the top of the Water Flosser handle. Press firmly until the tip clicks into place. To remove the tip from the handle, press the tip eject button and pull tip from the handle.

**CAUTION: Do not press tip eject button while unit is running.**
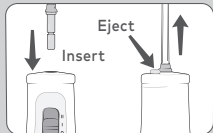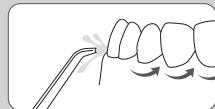
### Adjusting the Pressure Setting

You can adjust the pressure setting on your Water Flosser by simply sliding the pressure/power switch on the handle. Slide power switch up 1 click for Low (I) pressure and 2 clicks for High (II) pressure.

**Start with the lowest setting and increase to your liking.**

### Recommended Technique

Lean low over sink and place tip in mouth. Aim the tip toward gumline, and turn the unit ON (I).

For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth.


Continue until you have cleaned the inside and outside of both the upper and lower teeth.

Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. Keep unit upright during use for best results.

### When Finished

Turn the unit OFF (O). Empty any liquid left in the reservoir.

If desired, the reservoir may be removed from the power handle by simply pulling down on the reservoir to remove.

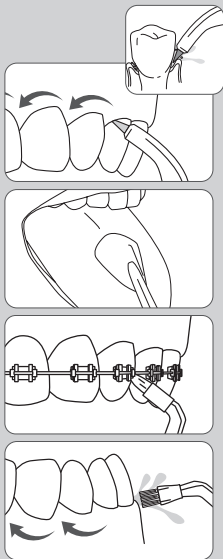**NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.**

## TIP USAGE

### Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.

To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket.

### Tongue Cleaner Tip

Set to lowest pressure setting, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure.

Increase pressure as you prefer.

### Orthodontic Tip

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

### Plaque Seeker™ Tip

To use the Plaque Seeker™ tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

### Using Mouthwash and Other Solutions

Your Waterpik™ Water Flosser can be used to deliver mouthwash and antibacterial solutions. For ease of adding other solutions remove the reservoir (see "When Finished" section) and add from top opening on reservoir. If you are adding additional water you may replace the reservoir and fill reservoir fully with the fill door. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

### Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior.

The water reservoir is removable for easy cleaning and is top rack dishwasher safe.

It is recommended that the reservoir is removed and cleaned weekly.

### Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 1 tablespoon of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

This process should be done every 1 to 3 months to ensure optimal performance.

### Battery Care

If the Water Flosser is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage.

Do not expose the battery or appliance to excessive temperatures.

### Battery Disposal/Recycling

Li-ion

This product contains a **rechargeable lithium ion battery**. The battery is not replaceable. At the end of product life, do not dispose of in household garbage or recycling bins. Please recycle at your nearest recycling center or take to hazardous waste collection point to comply with local and state regulations.

Dispose of used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire.

Please contact your local waste administrative department for any questions.

**Note: The rechargeable battery is not replaceable, remove only at end of product life. This process is not reversible and voids the product warranty.**

### Battery Removal

To remove the rechargeable battery for recycling, please discharge the battery fully before removing.



Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and surface on which you work.

1. Remove reservoir and wrap unit in towel. Hit with hammer until lower cover comes off.
2. Pull out battery, cut wires. Cut battery wires one at a time.
3. Recycle/dispose of the battery at your nearest recycling center or take to a hazardous waste collection facility to comply with local and state regulations.

### WARNING:

DO NOT cut both battery wires at the same time.

DO NOT disassemble, puncture, hit, squeeze or throw the battery into fire.
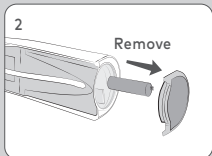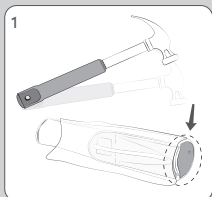
## Service Maintenance

Waterpik™ Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all your service needs go to [www.waterpik.ca](http://www.waterpik.ca). You may also call our toll-free customer service number at 1-888-226-3042. Refer to the serial and model numbers in all correspondence.

These numbers are located on the bottom of your unit.

Still have questions? Go to [www.waterpik.ca](http://www.waterpik.ca).

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion the product, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. This limited warranty excludes accessories or consumable parts such as tips, etc. Use with improper voltage will damage product and VOID warranty.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

WATER PIK, INC. DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

## **Limitation of Liability**

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence. This limited warranty is your sole remedy against Water Pik, Inc. and it shall not be liable for any damages, whether direct, indirect, incidental, special, consequential, exemplary or otherwise, including but not limited to property damage, arising out of any theory of recovery, including statutory, contract or tort.

## **Arbitration**

### BINDING ARBITRATION AND WAIVER OF RIGHT TO ASSERT MATTERS IN COURT AND RIGHT TO JURY TRIAL

By purchasing, installing, or using this product you agree that any controversy or claim arising out of or relating to your purchase, use, modification, installation, or removal of this product shall be settled by binding arbitration. The arbitrator's decision will be final and binding. Judgment on the award rendered by the arbitrator may be entered in any court having jurisdiction over the matter.

Further, you agree that you will not have the right to litigate any claim arising out of or in relation to your purchase, use, modification, installation, or removal of this product in any court or to have a jury trial on that claim.

Arbitration of any matters shall be administered by the American Arbitration Association and conducted by one arbitrator in accordance with the consumer

arbitration rules of the American Arbitration Association in effect at the time of filing of the demand. The venue for all such arbitrations shall be Colorado and any such arbitration shall be subject to the Federal Arbitration Act and Colorado law.

### WAIVER OF RIGHT TO PARTICIPATE IN CLASS ACTIONS OR PURSUE CLAIMS IN A REPRESENTATIVE CAPACITY

Further, by purchasing, installing, or using this product you agree that arbitration must be on an individual basis. This means neither you, nor we, may join or consolidate claims in arbitration by or against other purchasers or users of Water Pik, Inc. products or the general public. Additionally, you may not litigate in court or arbitrate any claims as a representative or member of a class or in a purported representative capacity on behalf of the general public, other purchasers or users of this product, other persons or entities similarly situated, or in a private attorney general capacity.

**IC Statement (Canada)**

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s).

Operation is subject to the following two conditions:

1. This device may not cause interference.
2. This device must accept any interference, including interference that may cause undesired operation of the device.

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**Water Pik, Inc., a subsidiary of /
une filiale de Church & Dwight Co., Inc.**